
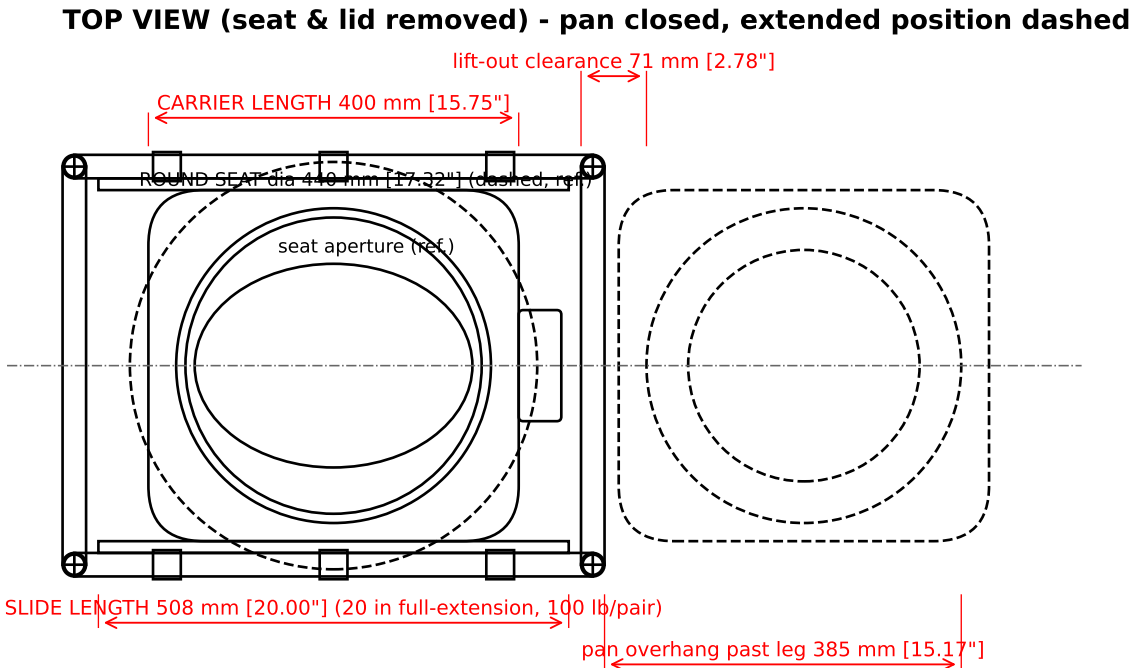


SHEET 1 - GENERAL ASSEMBLY: Folding commode with side slide-out pan (units mm [in], scale 1:1 in DXF)

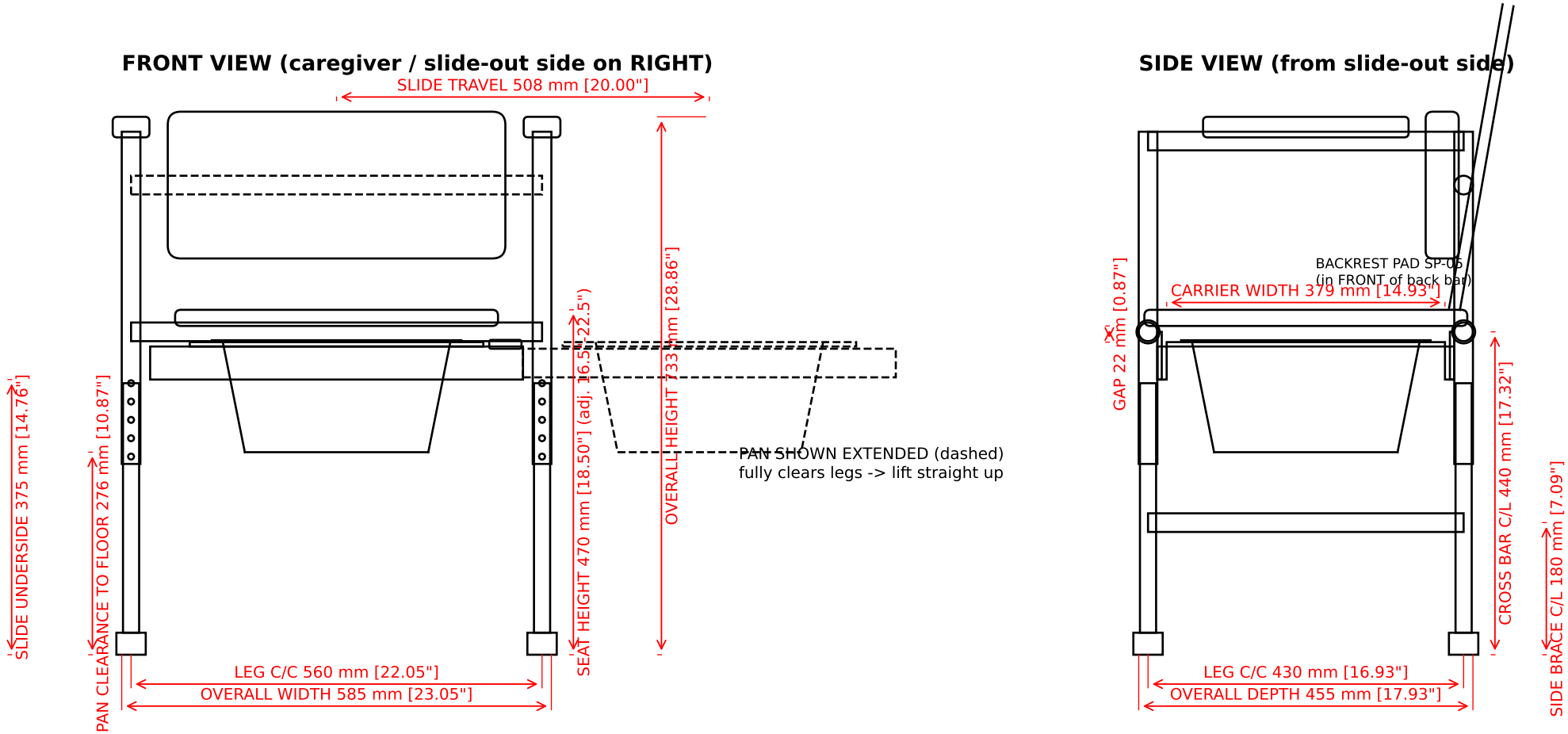

- NOTES
1. Base chair: standard folding aluminium commode, 1" (25.4 mm) x 1.5 mm wall tube, telescoping legs, 5-position height adjust.
 2. Slide-out pan module attaches to existing front & rear cross bars with 6x tube clamps (3 per bar) - no welding, retrofit-able.
 3. Two 20" full-extension ball-bearing drawer slides, 45 mm x 12.7 mm, stainless or zinc-plated, min. 100 lb (45 kg) rated per pair, with soft-close/detent.
 4. Carrier plate (blue) rides on slides; pan hangs by its rim in the carrier cutout and lifts straight out when extended.
 5. Pan slides out to the RIGHT of the seated user (caregiver side). Mirror assembly for left-hand version.
 6. Spring latch / magnetic catch holds carrier in CLOSED position; positive stop at full extension.
 7. All edges deburred. Static load: 300 lb (136 kg) user. Materials in Sheet 4 BOM.