




SHEET 1 - GENERAL ASSEMBLY: Folding commode with REAR slide-out pan (units mm [in], scale 1:1 in DXF)

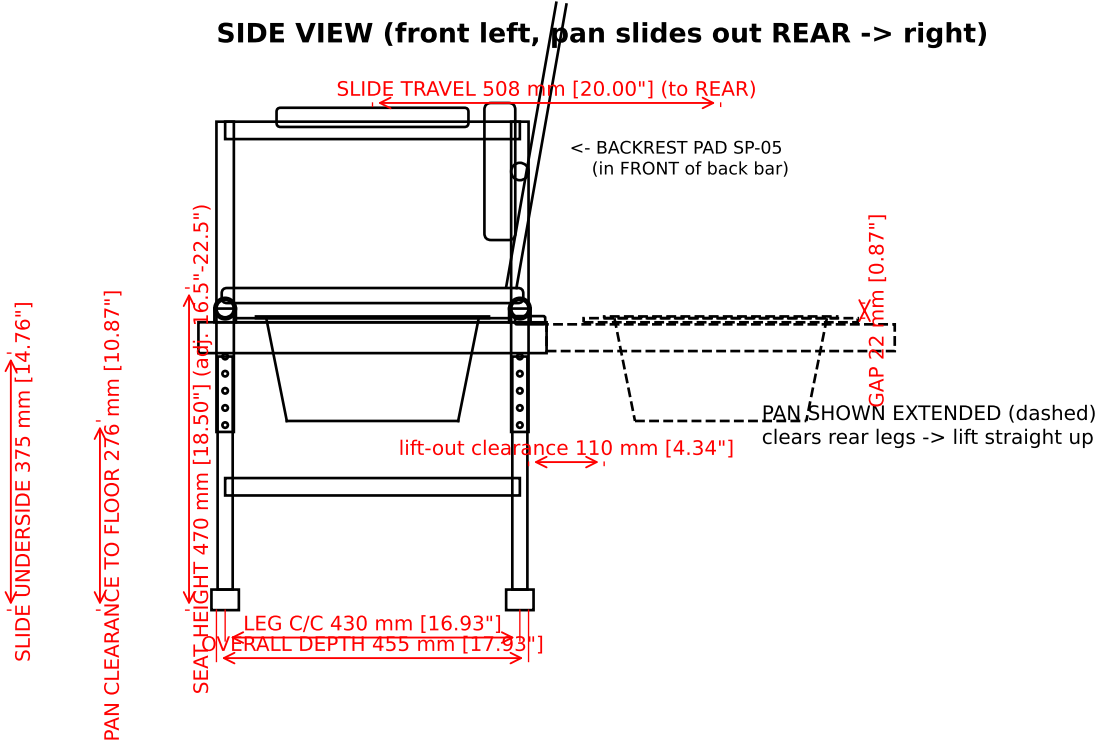

FRONT VIEW (slides on side rails, section through pan)

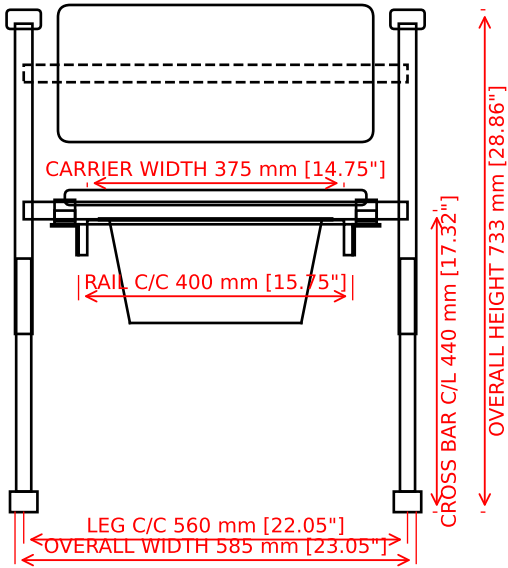

TOP VIEW (seat & lid removed) - pan closed, extended (rear) position dashed

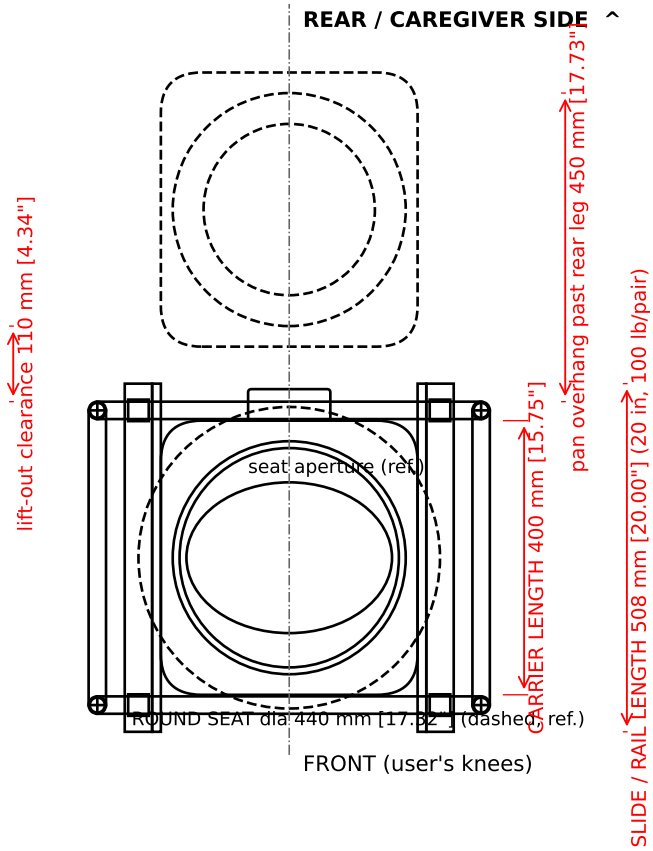

- NOTES
1. Base chair: standard folding aluminium commode, 1" (25.4 mm) x 1.5 mm wall tube, telescoping legs, 5-position height adjust.
 2. Slide-out pan module: two 40x40x3 aluminium angle rails (P/N SP-06) run front-to-back under the seat and clamp to the existing front & rear cross bars with 4x tube clamps (2 per rail) - no welding, retrofit-able.
 3. Two 20" full-extension ball-bearing drawer slides, 45 mm x 12.7 mm, stainless or zinc-plated, min. 100 lb (45 kg) rated per pair, with soft-close/detent, bolted to the rail webs.
 4. Carrier plate (blue) rides on slides; pan hangs by its rim in the round cutout and lifts straight out when extended.
 5. Pan slides out to the REAR of the seated user - caregiver stands behind the chair. Backrest pad (SP-05) and back bar are above the pan path; rear cross bar sits above the rails so nothing obstructs travel.
 6. Spring latch / magnetic catch holds carrier in CLOSED position; positive stop at full extension.
 7. All edges deburred. Static load: 300 lb (136 kg) user. Materials in Sheet 4 BOM.