


SHEET 4 - TUBE CLAMP BRACKET (P/N SP-04) & BILL OF MATERIALS (units mm [in])

CLAMP BRACKET - END VIEW
(2-piece saddle + foot, 30 mm long,
seen along cross bar)

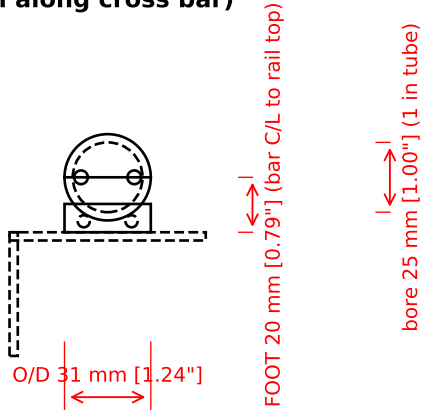

Material: 6061-T6 aluminium or glass-filled nylon.
2x M5 clamp screws (split line), 2x M4 screws down
into angle-rail top leg. Rubber liner 1 mm inside bore
to protect powder-coat & prevent slip.
Qty 4 (front & rear cross bar, at X = +/-220 mm).
RAIL SP-06: 40x40x3 6061-T6 angle, 508 mm long, web faces inboard;
slide cabinet member bolts to web (4x M4), qty 2.

BILL OF MATERIALS - slide-out pan module

ITEM	P/N	DESCRIPTION	MATERIAL	QTY
1	BASE	Folding commode chair, 1" alu frame, seat, lid, arms (existing / purchased)	Aluminium / PP	1
2	SLD-20	Full-extension ball-bearing drawer slide, 20" (508 mm), 45 mm, 100 lb/pair, detent	Zinc-plated or SS	2
3	SP-01	Carrier plate with formed flanges & cutout (Sheet 2)	6 mm HDPE or 2 mm 304 SS	1
4	SP-02	Removable pan, flanged rim, end grips (Sheet 3)	PP copolymer / 304 SS	1
5	SP-03	Snap-on splash lid for SP-02	PP copolymer	1
6	SP-04	Tube clamp bracket, 1" saddle with foot (this sheet)	6061-T6 or GF nylon	4
6a	SP-06	Angle rail 40x40x3 x 508 mm, drilled for slide + 2 clamps (this sheet)	6061-T6 aluminium	2
7	HDL-120	D-pull handle 120 mm, 12 mm dia	304 SS	1
8	LTCH	Magnetic catch 5 kg (or spring ball detent) + strike	SS / nylon	1
9	HW	M4x10 SS button-head screws + nyloc (slide to carrier flange)	A2 SS	12
10	HW	M5x20 SS screws + nyloc (clamp halves)	A2 SS	8
11	HW	M4x8 SS screws (clamp foot to rail; slide cabinet member to rail web)	A2 SS	16
12	STOP	Rubber bump stop 15 mm (full-extension stop on slide)	EPDM	2
13	SP-05	Backrest pad 460x200x45, wraps FRONT of back-bar, rear shell w/ 2 hook-&-loop straps	PU foam / vinyl, PP shell	1

- GENERAL
- Dimensions in mm [inches]. Untoleranced +/-1.0 mm. Break all edges 0.5 mm.
 - Frame contact parts must not damage powder-coat: use rubber-lined clamps.
 - Slide rails to be mounted parallel within 0.5 mm over 508 mm; check free travel with 20 kg in pan.
 - Design load: 136 kg (300 lb) seated user on base chair; 15 kg in pan on slides at full extension.
 - Finish: pan & carrier smooth, non-porous, cleanable with hospital disinfectants.
 - Pan extends to the REAR of the chair (caregiver behind user); no left/right handedness.